

Forced-Air Lab Oven


Scope of Application

The Forced-Air Lab Oven provides uniform thermal processing for laboratory drying, sterilization, and material storage. Optimized convection accelerates glassware/polymer drying while top exhaust ports managing moisture at $\leq 250^{\circ}\text{C}$.

Key Features


Simplified Operation

The LED display simplifies operation, making it easy to use and operate.

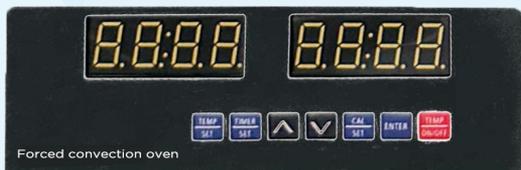

Uniform Distribution

The oven offers accurate temperature control and forced convection, ensuring uniform heat distribution throughout the chamber.

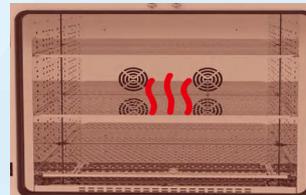

Venting Valve

The venting valve regulates the air exchange between the chamber and outside, ensuring optimal conditions for samples.



High Temperature Protection

The oven features high temperature protection with an alarm that sounds when the temperature exceeds the set point, ensuring safety for samples and personnel.



RS232 Interface

The RS232 interface allows for remote control by PC software, enabling convenient and efficient operation.





Specifications



| Model | MF-02G | MF-12G | MF-22G |
|------------------------------|-----------------------------------|-----------------|-----------------|
| Capacity (L) | 50 | 100 | 150 |
| Temp. Range (°C) | RT. + 10 ~ 250 | | |
| Temp.Stability at 100°C | 0.5 | 0.6 | 0.8 |
| Interior Dimension (W*D*H)mm | 370 x 355 x 420 | 480 x 460 x 500 | 500 x 480 x 625 |
| Exterior Dimension (W*D*H)mm | 573 x 622 x 771 | 600 x 690 x 810 | 660 x 750 x 920 |
| Number of Shelves | 2/8 | 2/12 | 2/14 |
| IP Code | IP20 | | |
| Net/ Gross Weight (kg) | 55/73.5 | 62/82 | 72.5/100 |
| Power Supply | 110V~/50/60Hz 220-240V~/50/60Hz | | |



| Model | MF-05G | MF-10G | MF-25G |
|------------------------------|-----------------------------------|-----------------|------------------|
| Capacity (L) | 50 | 100 | 250 |
| Temp. Range (°C) | RT. + 10 ~ 300 | | |
| Temp.Stability at 100°C | 0.5 | | |
| Interior Dimension (W*D*H)mm | 388 x 334 X 407 | 474 x 396 x 505 | 630 x 520 x 760 |
| Exterior Dimension (W*D*H)mm | 602 x 573 x 732 | 652 x 689 x 831 | 810 x 788 x 1031 |
| Number of Shelves | 2/4 | 2/12 | 2/14 |
| IP Code | IP20 | | |
| Net/Gross Weight (kg) | 55/73.5 | 62/82 | 78.5/106 |
| Power Supply | 110V~/50/60Hz 220-240V~/50/60Hz | | |

Optional(Foaming Glass Door)

